

Quick-RNA™ 96 Kit

RNA from any sample

Highlights

- 96-well spin-plate purification of total RNA (including small/microRNAs) from cells and tissue.
- DNA-free RNA is ready for Next-Gen Sequencing, RT/qPCR, etc. *DNase I is included.*

Catalog Numbers:

R1052, R1053



Scan with your smart-phone camera to
view the online protocol/video.



Table of Contents

Product Contents	01
Specifications.....	02
Product Description.....	03
Input Capacity & Average Yield Chart	04
Protocol	05
(I) Buffer Preparation	05
(II) Sample Preparation.....	06
Samples in DNA/RNA Shield, Cells, Tissue	06
(III) Total RNA Purification	07
Appendices	08
Sample Stabilization and Storage in DNA/RNA Shield	08
Samples in RNAprotect, RNAlater, UTM/VTM, etc.	08
Liquids/Reaction Clean-up	08
Homogenization with ZR BashingBead Lysis Tubes	09
Proteinase K Treatment	09
Protein Purification	09
Ordering Information	10
Complete Your Workflow.....	11
Troubleshooting Guide.....	12
Notes.....	13
Guarantee	17

Product Contents

Quick-RNA™ 96 Kit	R1052 (2 x 96 prep)	R1053 (4 x 96 prep)
RNA Lysis Buffer	100 ml (x2)	100 ml (x4)
RNA Prep Buffer	100 ml	100 ml (x2)
RNA Wash Buffer ¹ (concentrate)	48 ml (x2)	48 ml (x4)
DNase/RNase-Free Water	10 ml	30 ml
DNase I ² (lyophilized)	250 U (x4)	250 U (x8)
DNA Digestion Buffer	16 ml	16 ml (x2)
Silicon A™ Plate	2	4
Collection Plate	2	4
Elution Plate	2	4
Cover Foil	2	4
Instruction Manual	1	1

Storage Temperature - Store all kit components (i.e., buffers, columns) at room temperature.

Before use:

1 Add 192 ml 100% ethanol (208 ml 95% ethanol) to the 48 ml **RNA Wash Buffer** concentrate.

2 Reconstitute lyophilized **DNase I** with **DNase/RNase-Free Water**, mix by gentle inversion and store frozen aliquots:

#E1009-A (250 U), add 275 µl **water**

Specifications

- **Sample Sources** – Cells (animal, gram(-) bacteria), soft and easy-to-lyse tissue, samples in **DNA/RNA Shield™** or other preservation reagents, and enzymatic reactions (e.g., DNase I treated, Proteinase K treated). Not compatible with whole blood¹ and urine¹ samples.
- **Size** – Total RNA including small/microRNAs (≥ 17 nt).
- **Purity** – A_{260}/A_{280} & $A_{260}/A_{230} > 1.8$. RNA is ready for Next-Gen Sequencing, RT/qPCR, etc. Trace DNA can be removed by DNase I digestion (page 7).
- **Binding Capacity** – **Silicon A™ Plate** yield up to 10 μ g RNA/well.
- **Compatibility** – For samples stored in preservation reagents: **DNA/RNA Shield™**, RNaprotect®, Allprotect®, Universal transport medium/viral transport medium (UTM®/VTM®) and RNeasy Lysis Buffer™.
- **Elution Volume** – ≥ 25 μ l **DNase/RNase-Free Water**.
- **Equipment Needed (user provided)** – Centrifuge/rotor compatible with 96-well plates.


¹ For RNA purification from whole-blood and urine, see the Quick-RNA Miniprep Plus Kit (R1057, R1058) or the Quick-RNA MagBead Kit (R2132, R2133).

Product Description


The **Quick-RNA™ 96 Kit** provides a quick method for the isolation of high-quality total RNA (≥ 17 nt) from cells (animal, buccal, buffy coat, gram(-) bacteria) and soft, easy-to-lyse tissue.

The procedure combines a unique buffer system and spin-plate technology that results in high-quality total RNA (including small/microRNAs) that is ready for Next-Gen Sequencing, RT/qPCR, hybridization, etc.

High-Quality, *DNA-free* RNA

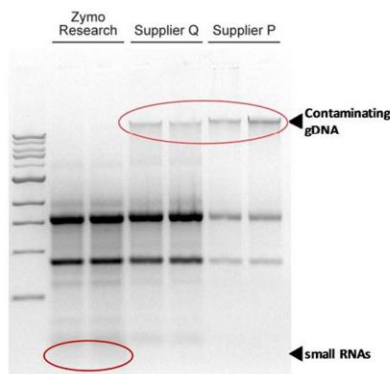

The **Quick-RNA™ Kits** yields high quality total RNA. High levels of genomic DNA contamination are present in the preps from Suppliers Q & P but not with the **Quick-RNA™ Kits**. Total RNA was isolated from human epithelial cells (sans DNase treatment).

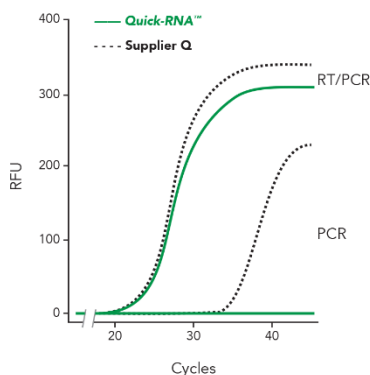

RNA isolated with the **Quick-RNA™ Kits** is DNA-free. Samples isolated with Supplier Q's kit are provided for comparison. Total RNA was isolated from 10^6 human epithelial cells (with in-column DNase treatments for both kits). Each amplification curve represents an average of three independent isolation experiments.

Input Capacity and Total RNA Yield

Input	Average RNA Yield	Kit Capacity
Cells	1 µg (per 10 ⁵ cells)	Up to 10 ⁶
HeLa	1.5 µg	
High Yield Tissue ^{1,2} (mouse)	≥ 3 µg (per 1 mg)	Up to 2 mg
Spleen	3-5 µg	
Liver	4-6 µg	
Low Yield Tissue ^{1,2} (mouse)	≤ 3 µg (per 1 mg)	Up to 5 mg
Brain, Heart	0.5-1.5 µg	
Muscle	0.5-2 µg	
Lung	1-2 µg	
Intestine	1-3 µg	
Kidney	2-3 µg	
Whole Blood ^{1,3}	(per 100 µl)	Up to 200 µl
Porcine	1-2 µg	
Human	0.2-1 µg	

1 Recommended: For protein-rich samples (e.g., tissue, blood cells, etc.) stored in DNA/RNA Shield™ (#R1100, sold separately), Proteinase K treatment can be performed using Proteinase K Set (#D3001-2-5, #D2001-2-20; sold separately) and PK Digestion Buffer (#R1200-1-5, #R1200-1-20; sold separately), see page 9.

2 Yield from tissue can vary due to other factors (i.e., organism type, physiological state, and growth conditions).

3 Yield from blood can vary based upon collection, sample preparation, donor, age, and/or health conditions.

Protocol

The protocol consists of: (I) Buffer Preparation, (II) Sample Preparation and (III) Total RNA Purification.

(I) Buffer Preparation

- ✓ Add 192 ml 100% ethanol (208 ml 95% ethanol) to the 48 ml **RNA Wash Buffer** concentrate.
- ✓ Reconstitute lyophilized **DNase I** with **DNase/RNase-Free Water**, mix by gentle inversion and store frozen aliquots: **#E1009-A (250 U)**, add 275 µl **water**

(II) Sample Preparation

- ✓ Perform all steps at room temperature and centrifugation at 10,000-16,000 x g for 30 seconds, unless specified.

Samples stabilized and stored in DNA/RNA Shield™ (cells, tissue, etc.)

1. If frozen, thaw homogenized sample in **DNA/RNA Shield™** to room temperature (20-30°C). Mix well by vortex.
2. Add an equal volume of **RNA Lysis Buffer** (1:1), mix well and proceed to purification, page 7.

Cells & Tissue (mammalian)

1. **Cells:**
 - a. If in suspension¹, pellet by centrifugation ($\leq 500 \times g$ for 1 minute), remove supernatant and resuspend cell pellet in **RNA Lysis Buffer** (see table below).
 - b. If adherent, remove liquid media from the culture container. Then add **RNA Lysis Buffer** directly to the monolayer (see table below). Remove cells from the culture surface by scraping, pipetting, etc.

Cells	Gram(-) Bacteria	Add RNA Lysis Buffer
$\leq 10^5$	-	$\geq 100 \mu\text{l}$
$\leq 10^6$	$\leq 10^8$	$\geq 300 \mu\text{l}$

To remove particulate debris, centrifuge and transfer the supernatant into a nuclease-free tube (not provided). Proceed to purification, page 7.

2. **Tissue**²:
 - a. Submerge an appropriate amount of fresh or frozen sample (see table below) into **RNA Lysis Buffer** and homogenize^{3,4}.

Tissue	Add RNA Lysis Buffer
High-yield ($\leq 2 \text{ mg}$) Low-yield ($\leq 5 \text{ mg}$)	$\leq 600 \mu\text{l}$

- b. To remove particulate debris, centrifuge and transfer the cleared supernatant into a nuclease-free tube (not provided). Proceed to purification, page 7.

1 If liquid/media cannot be removed, add ≥ 3 volumes **RNA Lysis Buffer** to 1 volume liquid sample (3:1) and mix well. Proceed to purification, page 7.
2 For examples of sample type input and average yield, see chart on page 4.
3 For efficient homogenization, bead beat samples with ZR BashingBead Lysis Tubes (S6012, S6003), sold separately. See Appendices (page 9) for bead beating parameters. Other types of homogenization can include mortar/pestle, dounce, syringe or tissue grinder, etc.
4 If no method of homogenization, alternatively tissue samples collected and stored in DNA/RNA Shield™ can be Proteinase K treated. See Appendices (page 9).

(III) Total RNA Purification

- ✓ Perform all steps at room temperature and centrifugation $\geq 2,500 \times g$ for 5 minutes (spin-plate).
- ✓ Do not use the **96-Well Plate Cover Foil** on the spin-plate during RNA Purification. If necessary, use an Air Permeable Sealing Cover (#C2011-8); sold separately.

1. Add 1 volume¹ ethanol (95-100%) to 1 volume sample lysed in **RNA Lysis Buffer** (1:1) and mix well.

Example: Add 300 μ l ethanol to 300 μ l mixture (sample lysed in **RNA Lysis Buffer**).

2. Transfer the mixture to a well of the **Silicon A™ Plate**² mounted on a **Collection Plate** and centrifuge³. Discard the flow-through.
3. **DNase I**⁴ treatment (recommended)
 - (D1) Add 400 μ l **RNA Wash Buffer** and centrifuge. Discard the flow-through.
 - (D2) For each sample/well to be treated, add 5 μ l **DNase I** (1 U/ μ l)*, 35 μ l **DNA Digestion Buffer** and mix by gentle inversion in an RNase-free tube (not included). Add 40 μ l directly to the column matrix of each well.
 - (D3) Incubate at room temperature (20-30°C) for 15 minutes. Proceed to step 4.
4. Add 400 μ l/well **RNA Prep Buffer** to the plate and centrifuge. Discard the flow-through.
5. Add 500 μ l/well **RNA Wash Buffer** to the plate and centrifuge. Discard the flow-through and repeat this step.
6. Mount the plate onto an **Elution Plate** and add 25 μ l **DNase/RNase-Free Water** directly to the matrix and centrifuge.

The eluted RNA⁵ can be used immediately or stored frozen.

Use the **96-Well Cover Foil** to prevent the eluate from evaporation.

1 To isolate large RNA species ≥ 200 nt, add 0.5 volume ethanol (95-100%) to 1 volume sample lysed in RNA Lysis Buffer (0.5:1) and mix well.

2 To process samples > 600 μ l, plate may be reloaded.

3 At this point, proteins can be purified from the flow-through (page 9).

4 Prior to use, reconstitute the lyophilized **DNase I** (Buffer Preparation, page 5). * Unit definition – one unit increases the absorbance of a high molecular weight DNA solution at a rate of 0.001 A₂₆₀ units/ml of reaction mixture at 25°C.

5 For complete removal of PCR (RT) inhibitors from plant, soil, and fecal samples, use the OneStep™ PCR Inhibitor Removal Kit (D6030).

Appendices

Sample stabilization and storage in DNA/RNA Shield™

Liquid samples (e.g., whole-blood): Add 3 volumes **DNA/RNA Shield™** (1X) to 1 volume sample (3:1). Mix well.

Solid samples (e.g., pelleted cells, tissue): Submerge sample (not to exceed 10% (v/v or w/v)) in **DNA/RNA Shield™** (1X) and homogenize (see Appendices, page 9).

Store samples in **DNA/RNA Shield™** at ambient temperature for ≥ 1 month or long term at frozen temperature. **DNA/RNA Shield™** is directly compatible with most guanidinium-based extraction methods (e.g., no need to remove reagent from the lysed/homogenized sample).

Samples in RNaProtect, Allprotect, RNaLater, UTM/VTM, saline or PBS

✓ RNaProtect®, Allprotect®: Add 3 volumes of **RNA Lysis Buffer** to 1 volume of liquid sample (3:1). Mix well and/or homogenize base on sample type (see Sample Preparation, page 6). Then, proceed to purification, page 7.

✓ RNaLater™:

- Cells - Pellet¹ by centrifugation at up to 5,000 x g and remove RNaLater (supernatant). Proceed to Sample Preparation, page 6.
- Tissue - Transfer into a new tube with forceps and remove any excess RNaLater™. Proceed to Sample Preparation, page 6.

Alternatively, for liquid samples from which RNaLater cannot be removed, add 1 volume of nuclease-free water (or PBS) to 1 volume liquid sample (1:1) and mix. Then add 4 volumes **RNA Lysis Buffer** to 1 volume sample/water (or PBS) mixture (4:1). Mix again and proceed to Total RNA Purification, page 7.

✓ Swab samples in UTM®/VTM®, saline or PBS: Remove swab and add 3 volumes of **RNA Lysis Buffer** to 1 volume sample (3:1). Mix well and proceed to purification, page 7.

Optional: To inactivate pathogens, store at room temperature prior to purification, add 1 volume **DNA/RNA Shield™** (2X concentrate) to 1 volume liquid sample (1:1) and mix well. Then proceed to Sample Preparation, Samples in **DNA/RNA Shield™**, page 6.

Liquids/Reaction Clean-up (DNase I treated RNA, in vitro transcriptions, etc.)

Add 150 μ l **RNA Lysis Buffer** to a ≥ 50 μ l liquid sample (3:1) and mix well. Proceed to purification, page 7, step 1.

¹ Different cells may react differently to centrifugation forces, and it is recommended to test the pelleting procedure with non-valuable samples first. Diluting RNaLater™ by 50% with cold PBS reduces solution density allowing for lower forces during cell pelleting (e.g., 500 x g).

(Appendices continued)

Homogenization with ZR BashingBead Lysis Tubes

- ✓ Recommended for complete and efficient homogenization of tough-to-lyse samples (e.g., tissue, plant, seed, microbes, etc.). Lysis tubes sold separately.
- ✓ For high-speed (e.g., MP Bio FastPrep-24, Bertin Precellys) and low-speed (e.g., Vortex Genie) homogenizers, bead-beating time optimization may be required.

Input	Tissue		Microbes
	Mammalian	Plant/Seed or Insect	Bacteria, Swab, Yeast, Stool/Soil
Cat. no. (lysis bead size)	S6003 (2.0 mm)	S6003 (2.0 mm)	S6012 (0.5 mm and 0.1 mm)
High-speed	30-60 sec	3-5 min	30-60 sec
Low-speed	3-5 min	15-20 min	5-10 min

Proteinase K Treatment

- ✓ Proteinase K treatment (optional) can be performed on protein-rich samples stored in **DNA/RNA Shield™** (#R1100, sold separately) (e.g., tissue, blood cells, etc.), using **Proteinase K Set** (#D3001-2-5, #D2001-2-20; sold separately) and **PK Digestion Buffer** (#R1200-1-5, #R1200-1-20; sold separately).
1. For every 300 µl of **DNA/RNA Shield™** sample, add 15 µl **Proteinase K** and 30 µl **PK Digestion Buffer**. Mix and incubate at room temperature (20-30°C) for ≥30 minutes (homogenized) or 2-5 hours (non-homogenized). Optimization may be required.
 2. To remove particulate debris from homogenized tissue, centrifuge and transfer the supernatant into a nuclease-free tube (not provided).
 3. Add **RNA Lysis Buffer** to the supernatant (1:1) and mix well. Proceed to purification, page 7.

Protein Purification: Acetone Precipitation of Proteins

- ✓ After the RNA binding (page 7, step 2), the protein content in the flow-through can be purified:
1. Add 4 volumes of cold acetone (-20°C) to flow-through (4:1) and mix.
 2. Incubate the samples for 30 minutes on ice.
 3. Centrifuge at max speed for 10 minutes. Discard the supernatant. Keep the pellet.
 4. Add 400 µl ethanol (95-100%) to the protein pellet. Centrifuge at max speed for 1 minute. Discard the supernatant.
 5. Air-dry the protein pellet for 10 minutes at room temperature.
 6. Resuspend and vortex the pellet in a buffer appropriate for downstream application (e.g., SDS-PAGE sample loading buffer).

1 Unit definition – one unit increases the absorbance of a high molecular weight DNA solution at a rate of 0.001 A₂₆₀ units/ml of reaction mixture at 25°C.

Ordering Information

Product Description	Catalog No.	Size
Quick-RNA™ 96 Kit	R1052	2 x 96 preps.
	R1053	4 x 96 preps.

Individual Kit Components	Catalog No.	Amount
RNA Lysis Buffer	R1060-1-50	50 ml
	R1060-1-100	100 ml
RNA Prep Buffer	R1060-2-25	25 ml
	R1060-2-100	100 ml
RNA Wash Buffer (concentrate)	R1003-3-24	24 ml
	R1003-3-48	48 ml
DNase/RNase-Free Water	W1001-10	10 ml
	W1001-30	30 ml
DNase I Set (lyophilized) (250 U supplied with DNA Digestion Buffer, 4 ml)	E1010	1 set
Silicon A™ Plates	C2001	2
Collection Plates	C2002	2
Elution Plates	C2003	2
96-Well Plate Cover Foil	C2007-2	2 pack
	C2007-4	4 pack

Complete Your Workflow

- ✓ For tough-to-lyse samples, use ZR BashingBead Lysis Tubes:

ZR BashingBead Lysis Tubes

2.0 mm beads #S6003	Plant/animal tissue
0.1 + 0.5 mm beads #S6012	Microbes
0.1 + 2.0 mm beads #S6014	Microbes in tissue/insects

- ✓ For isolation of DNA/RNA from any sample:

Quick-DNA/RNA Plus kits

Microprep Plus #D7005	From 1 cell and up
MagBeads #R2130	Automatable (Tecan, Hamilton, Kingfisher, etc.)

- ✓ For clean-up (purification) and concentration of any RNA sample. (e.g., from the aqueous phase of TRIzol® extractions) or from any enzymatic reaction (e.g., DNase I treated RNA):

RNA Clean & Concentrator kits


Microprep #R1013-R1014	DNase I Set included
MagBeads #R1082	Automatable (Tecan, Hamilton, Kingfisher, etc.)

- ✓ For NGS:

Zymo-Seq RiboFree Total RNA Library Prep kit

#R3000	12 preps
#R3003	96 preps

Troubleshooting Guide

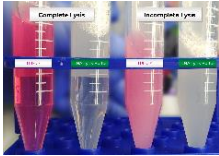Problem	Possible Causes and Suggested Solutions
Precipitation, viscous lysate	<p>Incomplete lysis and/or high-mass input:</p> <ul style="list-style-type: none"> - If precipitation occurs (upon adding ethanol to the lysate) or if the lysate is extremely viscous, increase the volume of DNA/RNA Shield™ and/or RNA Lysis Buffer to ensure complete lysis and homogenization until lysate is transparent (see image). 
Low purity (A_{260}/A_{230} nm, A_{260}/A_{280} nm)	<p>Sample handling:</p> <ul style="list-style-type: none"> - Ethanol and/or salt contamination. After centrifugation steps, carefully remove the column from the collection tube to prevent buffer carryover. Alternatively, blot emptied collection tube with a tissue or towel. - Make sure lysate and wash buffers have passed completely through the matrix of the column. This may require centrifuging at a higher speed and/or longer time. <p>Incomplete lysis and/or cellular debris:</p> <ul style="list-style-type: none"> - Increase the volume DNA/RNA Shield™ and/or RNA Lysis Buffer (proportionally) to ensure complete lysis and homogenization. Be sure to centrifuge and pellet any cellular debris then process the cleared lysate.
Low yield	<p>Sample input:</p> <ul style="list-style-type: none"> - Too much input or incomplete lysis/homogenization can cause cellular debris to clog or overload the column and result in compromised nucleic acid recovery. Use less input material and/or increase the volume DNA/RNA Shield™ and/or RNA Lysis Buffer. <p>High-protein content (blood, plasma/serum, etc.)</p> <ul style="list-style-type: none"> - Perform Proteinase K treatment to the sample prior to purification. See Appendices.
DNA contamination	<p>To remove DNA:</p> <ul style="list-style-type: none"> - Perform in-column DNase I treatment (page 9) or perform DNase I treatment post-purification, then re-purify the treated sample. - For future preps, increase the volume of DNA/RNA Shield™ and/or RNA Lysis Buffer to ensure complete lysis and homogenization of the sample.
RNA degradation	<p>To prevent RNA degradation:</p> <ul style="list-style-type: none"> - Immediately collect and lyse fresh sample into DNA/RNA Shield™ and/or RNA Lysis Buffer ensure stability. Homogenized samples can be stored frozen for later processing.

For technical assistance, please contact 1-888-882-9682 or email tech@zymoresearch.com

Notes

This image shows a single sheet of white paper with horizontal blue or grey ruling lines. The lines are evenly spaced and run across the width of the page. There are approximately 20 lines visible. The paper has a slight shadow on its right side, suggesting it's resting on a surface.

Notes

[illegible]

Notes

This image shows a single sheet of white paper with horizontal blue or grey ruling lines. The lines are evenly spaced and run across the width of the page. There are approximately 20 lines visible. The paper has a slight shadow on its right side, suggesting it's resting on a surface.

Notes

This image shows a single sheet of white paper with horizontal ruling lines. The lines are evenly spaced and run across the width of the page. There are no margins, text, or other markings on the paper.



100% satisfaction guarantee on all Zymo Research products, or your money back.

Zymo Research is committed to simplifying your research with quality products and services. If you are dissatisfied with this product for any reason, please call 1(888) 882-9682.

Integrity of kit components is guaranteed for up to one year from date of purchase. Reagents are routinely tested on a lot-to-lot basis to ensure they provide the highest performance and reliability.

This product is for research use only and should only be used by trained professionals. It is not for use in diagnostic procedures. Some reagents included with this kit are irritants. Wear protective gloves and eye protection. Follow the safety guidelines and rules enacted by your research institution or facility.

TM Trademarks of Zymo Research Corporation

Quick-RNA[®] is a registered trademark of Zymo Research Corporation. Other trademarks: TRI Reagent[®], TRIzol[®] (Molecular Research Center, Inc.), RNAp Protect[®], Allprotect[®], PAXgene[®] (Qiagen), RNAlater[®] (Ambion, Inc.), Universal (viral) transport medium (BD, Copan), Bioanalyzer (Agilent Technologies, Inc.).

*The **BEAUTY** of **SCIENCE** is to Make Things **SIMPLE**®*

