

RPA817Cp01 100µg

Recombinant Haptoglobin (Hpt)

Organism Species: *Capra hircus*; *Caprine* (Goat)

Instruction manual

FOR RESEARCH USE ONLY

NOT FOR USE IN CLINICAL DIAGNOSTIC PROCEDURES

12th Edition (Revised in Aug, 2016)

[PROPERTIES]

Source: Prokaryotic expression

Host: *E.coli*

Residues: Val19~Asn401

Tags: Two N-terminal Tags, His-tag and SUMO-tag

Subcellular Location: Secreted

Purity: > 90%

Traits: Freeze-dried powder

Buffer formulation: 20mM Tris, 150mM NaCl, pH8.0, containing 0.01% SKL, 5% Trehalose.

Original Concentration: 200µg/mL

Applications: Positive Control; Immunogen; SDS-PAGE; WB.

(May be suitable for use in other assays to be determined by the end user.)

Predicted isoelectric point: 8.0

Predicted Molecular Mass: 56.5kDa

Accurate Molecular Mass: 63kDa as determined by SDS-PAGE reducing conditions.

Phenomenon explanation:

The possible reasons that the actual band size differs from the predicted are as follows:

1. Splice variants: Alternative splicing may create different sized proteins from the same gene.
2. Relative charge: The composition of amino acids may affect the charge of the protein.
3. Post-translational modification: Phosphorylation, glycosylation, methylation etc.
4. Post-translation cleavage: Many proteins are synthesized as pro-proteins, and then cleaved to give the active form.
5. Polymerization of the target protein: Dimerization, multimerization etc.

[USAGE]

Reconstitute in 20mM Tris, 150mM NaCl (pH8.0) to a concentration of 0.1-1.0 mg/mL. Do not vortex.

[STORAGE AND STABILITY]

Storage: Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

Stability Test: The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

[SEQUENCE]

```

                VE TGSEAAAGSC PKAPEIANGH VEYSVRYQCD
KYYKLRAGNG VYTFNNKQWI NKDIGLQLPE CEEDVSCPEA PKIKNGYVEY
SVRYQCKTYY KLRTDGDGVY TFNSKKQWIN KNVGQQLPEC EAVCGKPKHP
VDQTQRIIGG SLDAGSFPW QAKMVSHHNL ISGATLINER WLLTTAKNLY
LGHTSDKKAK DITPTLRLYV GKNQLVEVEK VVLHPDHSKV DIGLIKREK
VPVNDKVMPI CLPSKDYVAV DRVGYSVGWG RNENFNFTGH LKYVMLPVAD
QDKCVKHYEG NNAPKNKTAT SPVGVQPILN ENTFCVGLSK YQEDTCYGDA
GSAFVVHDQE DDTWYAAGIL SFDKSCAVAE YGVYVKVTSI LDWVRKTIAN
N
    
```

[IDENTIFICATION]

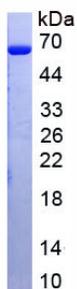


Figure. SDS-PAGE

[IMPORTANT NOTE]

The kit is designed for research use only, we will not be responsible for any issue if the kit was

used in clinical diagnostic or any other procedures.