


PRODUCT INFORMATION



BMS 191011

Item No. 22299

CAS Registry No.: 202821-81-6
Formal Name: 3-[(5-chloro-2-hydroxyphenyl)methyl]-5-[4-(trifluoromethyl)phenyl]-1,3,4-oxadiazol-2(3H)-one
MF: C₁₆H₁₀ClF₃N₂O₃
FW: 370.7
Purity: ≥95%
UV/Vis.: λ_{max}: 280 nm
Supplied as: A crystalline solid
Storage: -20°C
Stability: ≥2 years

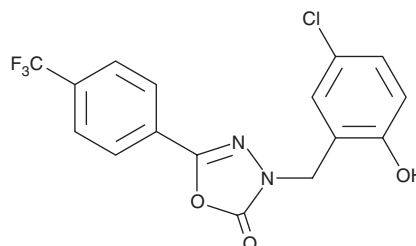

Information represents the product specifications. Batch specific analytical results are provided on each certificate of analysis.

Laboratory Procedures

BMS 191011 is supplied as a crystalline solid. A stock solution may be made by dissolving the BMS 191011 in the solvent of choice. BMS 191011 is soluble in organic solvents such as ethanol, DMSO, and dimethyl formamide (DMF), which should be purged with an inert gas. The solubility of BMS 191011 in ethanol is approximately 3 mg/ml and approximately 10 mg/ml in DMSO and DMF.

BMS 191011 is sparingly soluble in aqueous buffers. For maximum solubility in aqueous buffers, BMS 191011 should first be dissolved in DMSO and then diluted with the aqueous buffer of choice. BMS 191011 has a solubility of approximately 0.3 mg/ml in a 1:1 solution of DMSO:PBS (pH 7.2) using this method. We do not recommend storing the aqueous solution for more than one day.

Description

BMS 191011 is an activator of large-conductance calcium-activated potassium (K_{Ca}1.1/BK) channels that increases maximum potassium current to 126% of control in *X. laevis* oocytes expressing human K_{Ca}1.1 channels when used at a concentration of 1 μM.¹ *In vivo*, BMS 191011 (10-100 μg/kg, i.v.) increases the diameter of retinal arterioles without affecting blood pressure or heart rate in rats, an effect that is reversed by the K_{Ca}1.1 channel blocker iberiotoxin (Item No. 14608).²

References

1. Romine, J.L., Martin, S.W., Meanwell, N.A., et al. 3-[(5-Chloro-2-hydroxyphenyl)methyl]-5-[4-(trifluoromethyl)phenyl]-1,3,4-oxadiazol-2(3H)-one, BMS-191011: Opener of large-conductance Ca²⁺-activated potassium (maxi-K) channels, identification, solubility, and SAR. *J. Med. Chem.* **50**(3), 528-542 (2007).
2. Mori, A., Suzuki, S., Sakamoto, K., et al. BMS-191011, an opener of large-conductance Ca²⁺-activated potassium channels, dilates rat retinal arterioles *in vivo*. *Biol. Pharm. Bull.* **34**(1), 150-152 (2011).

WARNING

THIS PRODUCT IS FOR RESEARCH ONLY - NOT FOR HUMAN OR VETERINARY DIAGNOSTIC OR THERAPEUTIC USE.

SAFETY DATA

This material should be considered hazardous until further information becomes available. Do not ingest, inhale, get in eyes, on skin, or on clothing. Wash thoroughly after handling. Before use, the user must review the [complete](#) Safety Data Sheet, which has been sent via email to your institution.

WARRANTY AND LIMITATION OF REMEDY

Buyer agrees to purchase the material subject to Cayman's Terms and Conditions. Complete Terms and Conditions including Warranty and Limitation of Liability information can be found on our website.

Copyright Cayman Chemical Company, 11/12/2021

CAYMAN CHEMICAL

1180 EAST ELLSWORTH RD
ANN ARBOR, MI 48108 · USA

PHONE: [800] 364-9897
[734] 971-3335

FAX: [734] 971-3640

CUSTSERV@CAYMANCHEM.COM
WWW.CAYMANCHEM.COM